

MICRONIC

MICRONIC | SAMPLE STORAGE | TUBES

0.75ML V-BOTTOM TUBE RANGE **96-WELL CONFIGURATION**



DO YOU NEED MORE SPACE FOR YOUR STORAGE AND STILL BE ABLE TO USE SINGLE TUBES IN A 96-WELL CONFIGURATION? IF SO, THEN MICRONIC OFFERS YOU THE POSSIBILITY TO SAVE BETWEEN 35% AND 50% OF YOUR SPACE. WHETHER YOU CHOOSE THE 0.75ML NON-CODED RANGE OR CODED RANGE, THE 0.75ML TUBE DESIGN WILL BE IDENTICAL AND WILL CONTAIN THE SAME FEATURES.

Micronic offers a comprehensive range of non-coded, alphanumeric coded and 2D coded (Data-Matrix and TraXis) tubes. Manufactured to the industry leading strict tolerances, the precision and tube-to-tube consistency of Micronic products maximises the operational uptime when used with automated cherry picking and handling systems. The tubes are manufactured of the highest-purity polypropylene, which resists many organic solvents.

The 2D codes are encrypted with a laser in the polypropylene code area, in turn guaranteeing the highest resistance to solvent attack and ensuring reliable code reading.

The tubes have an inner V-bottom (also available with U-bottom). The gross volume of the tube is 0.75ml, with net working volume ranges from 0.44ml to 0.52ml, depending on the capping style.

The 0.75ml tubes are compatible with the Micronic Plug Style Caps or Screw Caps (TPE or EVA Caps or Screw Caps also fit Screw Cap tubes).

The 0.75ml tubes are available in bulk, refill system-96 and Lobarack-96, which features 'Twist-Lock' to prevent turning during Screw Capping.

Available features

- Amber tubes for light sensitive samples
- Sterile, gamma radiated
- Snap Tubes, locks tubes into rack wells to prevent sample loss from overturned racks
- Precapped, with TPE Caps or Screw Caps

For more information:
e-mail to sales@micronic.com


Micronic Europe B.V.

(Asia, Africa, Europe and Oceania)
Phone: +31 (0)320 277 070

Micronic America LLC


(North, Central and South America)
Phone: +1 484 480 3372

Please visit our website: www.Micronic.com

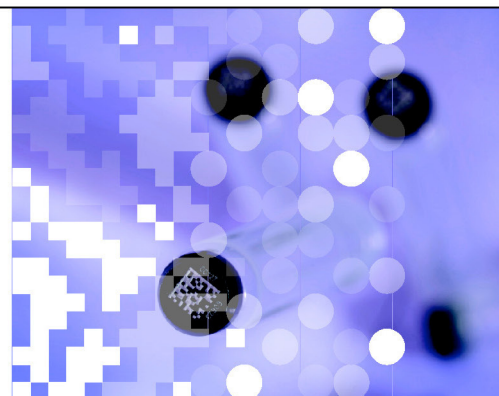


INFORMATION

- Gross volume 0.75ml
- Net volume 0.52ml when used with a Plug Style cap
- Net volume 0.44ml when used with a Screw Cap
- Available with an inner V-bottom
- Space saving shorter tube
- Automation friendly
- 2D codes are permanently laser-etched in the polypropylene
- 2D codes also have a human-readable code
- Resists high chemical solvents
- Class 7 clean-room production
- RNase/DNase and Pyrogene free
- 100% quality control on all 2D codes

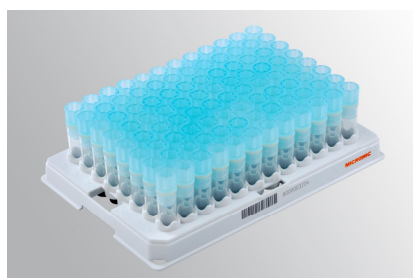

MICRONIC

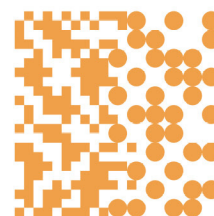


MICRONIC

MICRONIC | SAMPLE STORAGE | TUBES

SPECIFICATIONS

MATERIAL	HIGHEST PURITY PP
GROSS VOLUME	0.75ML
NET VOLUME	0.44ML SCREW CAP/0.52ML PLUG STYLE CAP
RESISTANCE	HIGH CHEMICAL SOLVENTS
CODING	NON-CODED, ALPHANUMERIC PRINTED, OR 2D CODED (DATA-MATRIX OR TRA Xis)

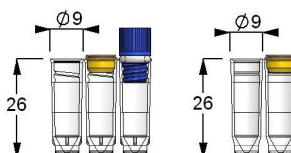

ORDERING INFORMATION


See also the leaflet "Line of Screw Cap Storage tubes"

0.75ml Non-coded Screw Cap tubes V-bottom without cap		
MP32069	0.75ml NC SC tubes V in bulk	box/960
MPW32070	0.75ml NC SC tubes V in Loborack-96	box/10
MPW32071	0.75ml NC SC tubes V in Loborack-96 Low Cover for TPE C	box/10
0.75ml Alphanumeric coded tubes V-bottom without cap		
MP42067	0.75ml AN SC tubes V in refill system-96	box/10
MPW42068	0.75ml AN SC tubes V in Loborack-96	box/10
0.75ml 2D Data-Matrix coded Screw Cap tubes V-bottom without cap		
MP52349	0.75ml 2D Data-Matrix C SC tubes V in bulk	box/960
MP52350	0.75ml 2D Data-Matrix C SC tubes V in refill system-96	box/10
MPW52351	0.75ml 2D Data-Matrix C SC tubes V in Loborack-96	box/10
0.75ml 2D TraXis coded Screw Cap tubes V-bottom without cap		
MP52639	0.75ml 2D TraXis C SC tubes V in bulk	box/960
MP52640	0.75ml 2D TraXis C SC tubes V in refill system-96	box/10
MPW52638	0.75ml 2D TraXis C SC tubes V in Loborack-96	box/10
0.75ml SC tubes V-bottom pre-capped with Screw Cap Transparent Blue		
MP32069-Z6	0.75ml NC SC tubes V with SC in bulk	box/960
MPW32071-Z6	0.75ml NC SC tubes V with SC in Loborack-96	box/10
MPW42068-Z6	0.75ml AN SC tubes V with SC in Loborack-96	box/10
MP52349-Z6	0.75ml 2D Data-Matrix C SC tubes V with SC in bulk	box/960
MPW52351-Z6	0.75ml 2D Data-Matrix C SC tubes V with SC in Loborack-96	box/10
0.75ml SC Tubes V-bottom pre-capped with Pierceable TPE Cap Natural		
MP32008-X01	0.75ml NC SC tubes V with TPE C in Loborack-96	box/10
MPW52020-X01	0.75ml 2D Data-Matrix SC tubes V with TPE C in Loborack-96	box/10

INFORMATION

- NC = Non-coded
- AN = Alphanumeric coded
- Data-Matrix C = 2D Data-Matrix coded
- TraXis C = 2D Traxis coded
- SC = Screw Cap
- V = V-bottom
- TPE = Thermo Plastic Elastomer
- TPE C = Thermo Plastic Elastomer Cap

Micronic has a broader product range than the products listed under ordering information. For our complete product portfolio, please send your e-mail to sales@micronic.com or contact us by phone.

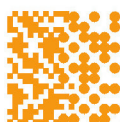

MICRONIC